

URETHRAL CATHETERISATION DEVICE (UCD®)

Video Demonstation & Feedback

We'd love to hear from you and support product adoption at your institution. Simply scan the QR Code to watch a product demonstration and access our feedback form, or visit <https://scnv.io/dVOC>



Device Description

The innovative Urethrotech UCD® (Urethral Catheterisation Device) is composed of a 3-way Silicone Foley Catheter with integrated introducer guidewire for safe bladder access to prevent Catheterisation-associated Urethral Injury (CAUI). The integrated removable guidewire is made of Nitinol with hydrophilic coating to facilitate the atraumatic introduction of the catheter through the urethra into the bladder to negotiate any prostatic curve or tight external urethral sphincter. Also suitable for irrigation of the bladder after guidewire removal.

- >> Completely 100% Silicone 3-way Foley catheter
- >> Clear for easy visualisation of urine flow or debris
- >> Radiopaque blue tip and stripe for easy visualisation during scans
- >> Two large opposing eyes for efficient urine drainage
- >> Soft atraumatic integrated hydrophilic Nitinol guidewire
- >> Guidewire exits at tip of catheter for ideal catheter guidance
- >> Secure guidewire Luer-lock for lubrication of guidewire
- >> Attached plug ready to seal guidewire channel after guidewire removal
- >> Symmetrical balloon shape for added reliability
- >> Latex Free

Urethral Catheterisation Device (UCD®)

3-Way 100% Silicone Foley Catheter with integrated guidewire

PACKAGING						1box/each
Unit of	Net	Volume	Shipping Dimensions			
Measure	Conversion	Weight (g)	(cubic m³)	L x W x H (cm)	GTIN-13	
Box (BX)	1.0 Each	111	0.00072	52.2 x 9.2 x 1.5	5060639870004	
Each (EA)	1.0 Each	111	0.00072	52.2 x 9.2 x 1.5	5060639870004	

Size

Item REF No.	Size	Balloon	Length	NPC code	NHS/CCSD Codes	Pkg.
UCD-0035-090-0014	14 Fr	10 mL	90cm	FST.....	M478,M371/M4820	1/bx

Product Specifications

HPIS CODE	730.10.40.20
BALLOON SIZE (CC)	10
CATHETER MATERIAL	Silicone, 100% Silicone
COATING	No
FRENCH SIZE (FR)	14
LATEX FREE	Yes
PRODUCT STYLE	3-Way Foley Catheter, UCD
INTEGRATED GUIDEWIRE (CM)	90
SILVER	No
TEMPERATURE SENSING	No
UNSPSC	42142702

Further Information: www.urethrotech.com

Product Training: training@urethrotech.com




Customer Service: sales@urethrotech.com

Direct Order: sales@urethrotech.com

Urethrotech Ltd
PO Box 1605, KINGSTON UPON THAMES KT1 9FW
Tel: 0800 051 7219 <https://urethrotech.com>

DISCLAIMER: Private Practice CCSD codes and NHS National Tariff Procedure and HRG-codes provided by Urethrotech are intended as general guidelines only. Urethrotech does not guarantee coverage or reimbursement of any products. The customer must address all coverage and reimbursement issues (including correctness and accuracy of codes) with their individual payers. It is the customer's responsibility to ensure the accuracy and appropriateness of each claim the customer submits, in accordance with all applicable payer requirements. We reserve the right to correct any errors that may occur within this brochure. © 2025 Urethrotech, Ltd. Urethrotech and UCD are registered trademarks of Urethrotech Ltd.

Ref: UCDCAT.UK.06/25

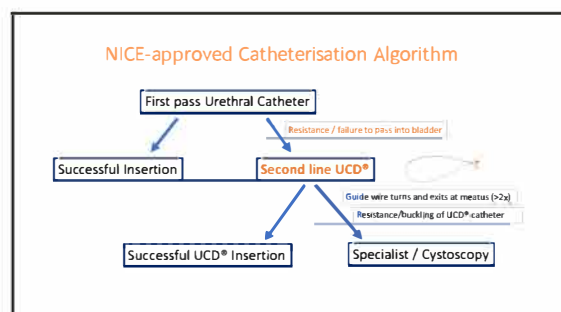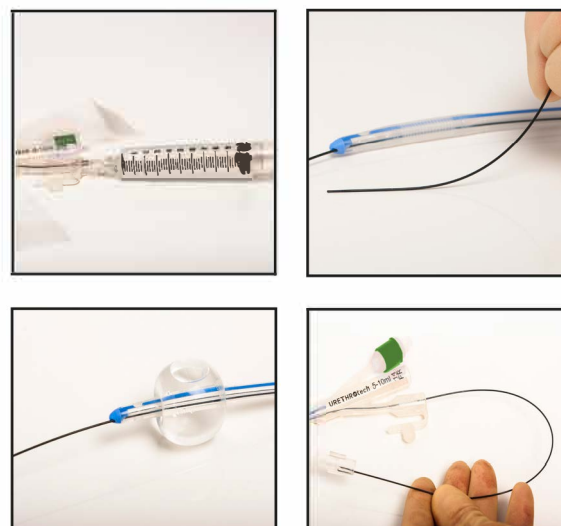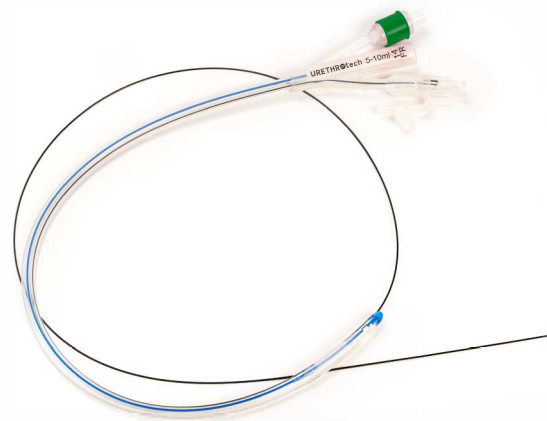

Order via NHS Supply Chain:

